


REPORT OF ANALYSIS

For: (35432) OPEN FARM - Dog-OF Pork Kibble
257221004

Analysis	Level Found		Reporting			Analyst- Date	Verified- Date
	As Received	Units	Limit	Method			
Sample ID: 257221004-Micro		Lab Number: 13899209		Date Sampled: 2022-09-15			
Listeria	negative	org/125g	1	RapidChek/AOAC RI 020401	kkb0-2022/09/18	jzh4-2022/09/18	
E. coli (generic)	n.d.	cfu/g	10	AOAC 2018.13	Kat6-2022/09/17	snl7-2022/09/17	
Salmonella	negative	org/375g	1	AOAC 2003.09; AFNOR QUA 18/03-11/02	snl7-2022/09/17	jzh4-2022/09/19	

All results are reported on an AS RECEIVED basis., n.d. = not detected , cfu = colony forming unit

For questions please contact:

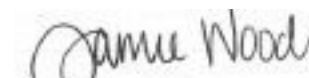

Jamie Wood
Account Manager
jwood@midwestlabs.com (402)590-2964

The result(s) issued on this report only reflect the analysis of the sample(s) submitted.

Our reports and letters are for the exclusive and confidential use of our clients and may not be reproduced in whole or in part, nor may any reference be made to the work, the results, or the company in any advertising, news release, or other public announcements without obtaining our prior written authorization.

REPORT OF ANALYSISFor: (35432) OPEN FARM - Dog-OF Pork Kibble
257221004**Detailed Method Description(s)****E. coli and Total Coliform using 3M Petrifilm**

Sample analysis follows MWL MI 292 which is based on AOAC 2018.13. A representative sample is obtained and added to phosphate buffer. Aliquots of the sample are withdrawn and placed on Petrifilm plates. The plates are incubated for 18 to 24 hours. After incubation, the plates are counted to determine the number of generic E. coli and total coliforms present. The color of the colony and the presence of gas differentiate a generic coliform from E. coli. The levels are reported as colony forming units (cfu).

Listeria Lateral Flow

Samples are analyzed following MWL MI 194 which is based on the RapidChek Listeria User Guide. A representative sample is obtained and combined with a selective growth media. It is incubated for 40-48 hours. After incubation, an aliquot is heated for 10 minutes, and a test strip for Listeria detection is used. Results are reported as negative or presumptive positive. This procedure does not speciate Listeria.

Salmonella PCR

Sample analysis follows MWL MI 180 which is based on AOAC 2003.09. A representative sample is obtained and combined with Buffered Peptone Water. The sample is incubated for 16 hours. An aliquot of enriched sample is transferred to BHI and incubated for three hours. The enriched media is then analyzed by PCR for Salmonella detection. Results are reported as negative or presumptive positive.

The result(s) issued on this report only reflect the analysis of the sample(s) submitted.

Our reports and letters are for the exclusive and confidential use of our clients and may not be reproduced in whole or in part, nor may any reference be made to the work, the results, or the company in any advertising, news release, or other public announcements without obtaining our prior written authorization.